

**Cetis**  | **APPLIANCES**  
Teledex • TeleMatrix • Scitec

# **User Guide**

## **Weight Scale (S-100)**


Thank you for purchasing the Cetus body weight scale. Please read this user guide before using the scale. If you have any questions, please contact us at [www.cetisappliances.com](http://www.cetisappliances.com)

## PACKING LIST

The packaging box includes the following items:

- Weight Scale (**batteries not included**)

### Part Names:

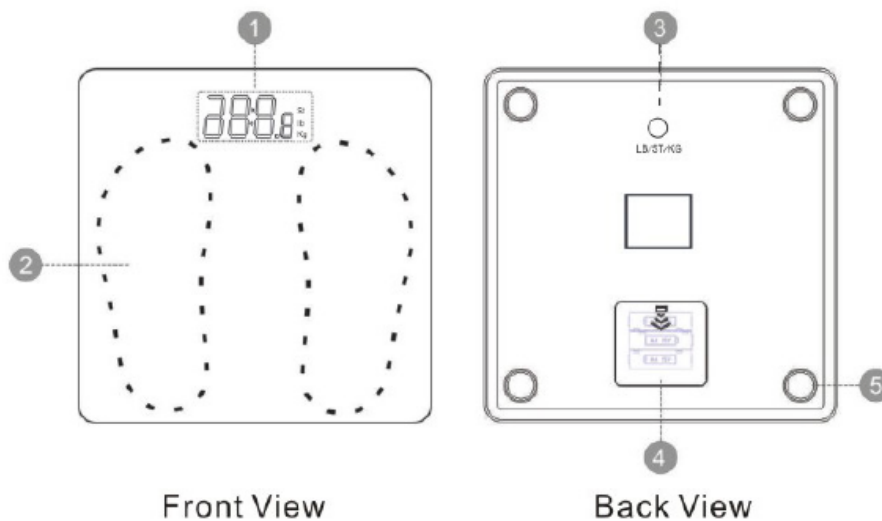

1. Display area
2. Standing area
3. Button to switch weight units
4. Battery compartment
5. Non-slip feet

Note: The button to switch weight units cycles through lb/St/Kg.

## PRODUCT FEATURES

- Patented series design
- 8-mm ultra-white tempered glass, compression tested 1250 times
- Ultra-thin anti-skid design to prevent rollovers
- White light number display for use at night
- Number display has delay for longer reading time
- One button for switching weight units
- Available in black or white



## **OPERATING INSTRUCTIONS**

1. Remove the battery compartment cover, install 3 AA alkaline batteries.
2. Using the weight unit switching button, select the appropriate weight unit. To switch weight units, gently press the front side to wake up the display area and then use the weight unit switch button.
3. Place the body weight scale on the floor, number display facing up.
4. Step on top of the body weight scale, so that display numbers are visible above your toes. Make sure both feet are on the scale. Wait until the display numbers stop shifting. Once the display numbers no longer change, this is your weight.
5. Step off the scale. The display area automatically turns off in 10 seconds and enters the standby mode. To wake it up, press the scale surface gently.

## Precautions

1. Do not stack objects on the body weight scale when it is not in use.
2. To prevent slipping, stand on top of the scale with both feet.
3. To ensure an accurate reading, stand on the scale without moving or shifting your weight.
4. To prevent slipping, do not stand on the scale if it is wet or if your feet are wet.
5. Do not hit or kick the scale. Do not drop or throw the scale.
6. Use the body weight scale in a dry area, preferably one with low humidity.

### COMMON FAILURES AND SOLUTIONS

Failure	Solution
There is no display.	Check if the batteries are correctly installed and if the batteries are rusty. If they are normal and correctly installed then notify the maintenance department about the issue.
Battery symbol is displayed.	Replace the batteries.

## TECHNICAL PARAMETERS

- Product name: Body Weight Scale
- Product model: S-100
- Rated power supply: 3 AA (#5) batteries **(batteries not included)**
- Low voltage indication: When the voltage is less than 3.5V, the battery symbol "Lo" is displayed when the scale is on
- Weighing range: 3.0–180 kg
- Display: LED SMD white light (weighing value, weight unit, battery charge level)
- Dimensions: 292 × 292 × 30 mm, Weight: 2350 g
- Certification: ETL & CE certification for safety
- Available in Black or White



The background of the entire page is an abstract, artistic composition of overlapping, translucent blue waves and ribbons. These shapes flow and swirl across the white background, creating a sense of movement and depth. The colors range from light, airy blues to deeper, more saturated tones, with some areas appearing as thin, delicate lines and others as thicker, more substantial folds.

**[www.cetisappliances.com](http://www.cetisappliances.com)**

**4975 N. 30th Street  
Colorado Springs, CO 80919**

**info\_USCA@cetis.com  
Telephone +1 719 638 8821  
Fax +1 719 638 8757**

© 2022 Cetis, Inc. Cetis, Inc. is a U.S. corporation, incorporated in the State of Delaware and headquartered in Colorado. Product specifications and descriptions in this document subject to change without notice. Cetis® is trademarked or registered trademarks of Cetis, Inc.

**CTS-APP-UG-WS-112022**